


Your path to predictable

Everything you need to take your practice to the next level your way

Unparalleled technology is here to help you grow your practice efficiency and patient satisfaction:

- Faster, smoother scanning with superior visual output for enhanced metal and edentulous scanning
- Easier handling and advanced ergonomics

System includes:
Scanner, Three year warranty, and Laptop availability

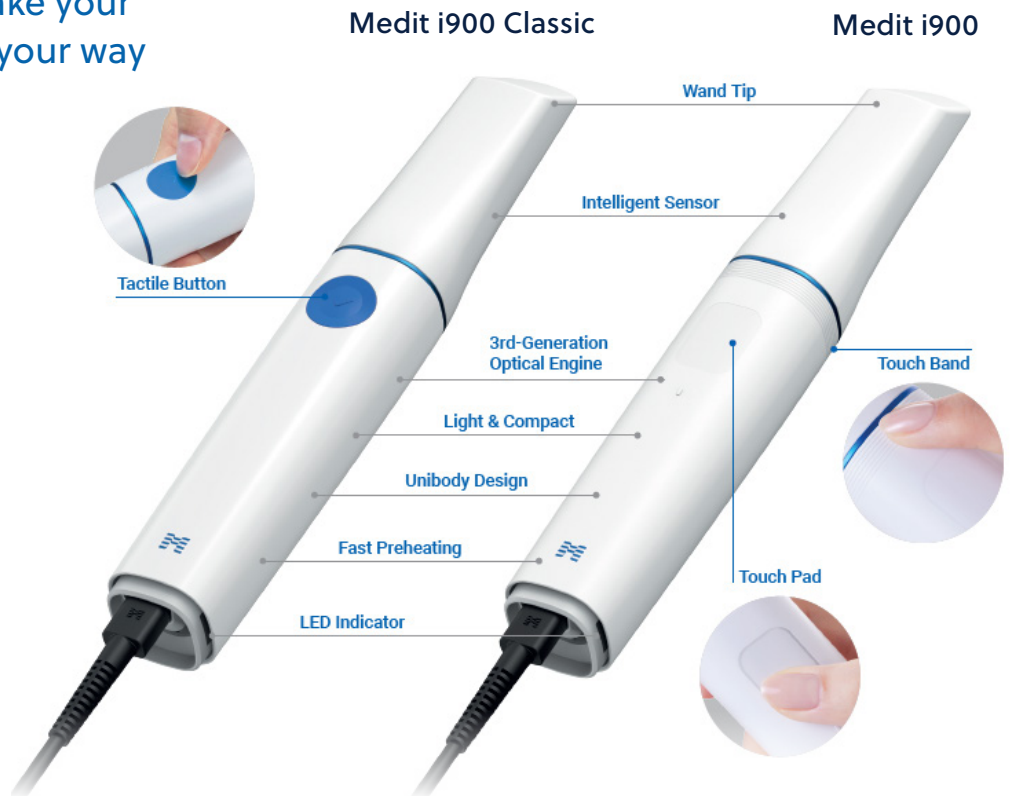

Medit Link

Secure cloud-based 3D data sharing platform

- Facilitates communication
- Encourages collaboration with dental partners
- Seamless Integration with **RealGUIDE™** and numerous other systems.

Medit Apps

A wide range of features to support restorative outcomes:

- Scan and color filtering
- Margin line creation
- Ortho simulation
- Temporaries and crown fit
- A.I. scan body matching



Unleash your clinic's full potential

Medit Intraoral Scanner Lineup Comparison



Find Your Perfect Scanner

Category		i900 classic	i900	i700 wireless	i700	i600
Model						
Scanning Technology	Imaging Technology	3D-in-motion video technology / 3D full color streaming capture				
	Anti Fogging Technology	Adaptive anti-fogging				
Accuracy	Full Arch	10.9 µm ± 0.98*				
Handpiece	Interface Type	Button	Touch	Button		
	Dimensions	223 x 37 x 35 mm		313 x 44 x 47 mm	248 x 44 x 47 mm	
	Total Weight	165g		328g	245g	
Tip	Scan Area	18 x 15 mm		15 x 13 mm		
	Autoclavable	Sterilize for 10 minutes at 135°C (275°F) at gravity type and dry for 30 minutes. Sterilize for 4 minutes at 134°C (273°F) at pre-vacuum type and dry for 20 minutes.				
	Reversible Tip	No		Yes		Yes
Special Mode	Remote Control Mode	Yes	Yes			No
	UV-LED Disinfection	Yes			No	
	Smart Shade Guide	Yes			No	
	HD Mode	Yes			No	
MAC OS Compatibility		Yes			No	
Connectivity	Connection Type	Wired USB 3.0				
	Plug and Scan	Yes				
	Detachable Cable	Yes				

* May vary depending on scan object

For more information, visit ZimVie.com

ZimVie
4555 Riverside Drive
Palm Beach Gardens, FL 33410
1-800-342-5454
Phone: +1-561-776-6700
Fax: +1-561-776-1272



ZimVie Education



Medit User Group



Medit Academy



Medit Help Center



Unless otherwise indicated, as referenced herein, all trademarks are the property of ZimVie; or and affiliate; and all products are manufactured by one or more of the dental subsidiaries of ZimVie Inc. (Biomet 3i, LLC, Zimmer Dental, Inc., etc.) and marketed and distributed by ZimVie Dental and its authorized marketing partners. Intraoral scanners and accessories are manufactured by Medit Corp. For additional product information, please refer to the individual product labeling or instructions for use. Product clearance and availability may be limited to certain countries/regions. Medical device class I (EU). Medit scanners are intraoral 3D scanners intended to record topographical characteristics of teeth and surrounding tissues digitally. These medical devices can not be reimbursed by the French public insurance. This material is intended for clinicians only and does not comprise medical advice or recommendations. Distribution to any other recipient is prohibited. This material may not be copied or reprinted without the express written consent of ZimVie.
ZV2397 GLOB REV B 05/25 ©2025 ZimVie. All rights reserved.

